

Configuring the MF2 for Closed Field Operation

Important: Closed Field Adapter, and CF Line Filter used only for Closed Field Configuration!


1

Front of Speaker



Holes for Hex Screws

Back of CF Adapter

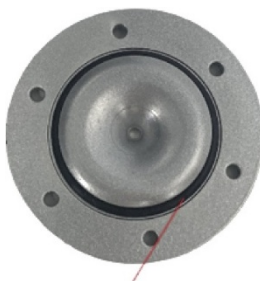

Black O-ring

Ensure Black O-ring is in place on back of CF Adapter, as shown.

Attach the CF Adapter to the front of the speaker using six of the provided 1/4 x 4-40 hex screws.

2

Speaker with CF Adapter



Ensure the black O-ring is in place on the back of the CF adapter as shown above.

Gently insert the tube into the narrow opening on the face of the CF Adapter.

3

CF Line Filter



For Closed Field Configuration Only

Attach the CF filter to the provided RCA cable.

4



Back of Speaker
(shown with included stand)


Amplifier Line In

Attach a BNC to RCA adapter to the BNC amplifier port of your source device.

Connect the MF2 to the amplifier using the RCA cable [with CF filter attached].

See the System 3 Manual for more information on the MF2 Speaker.

MF2 Closed Field Speaker Assembly Diagram

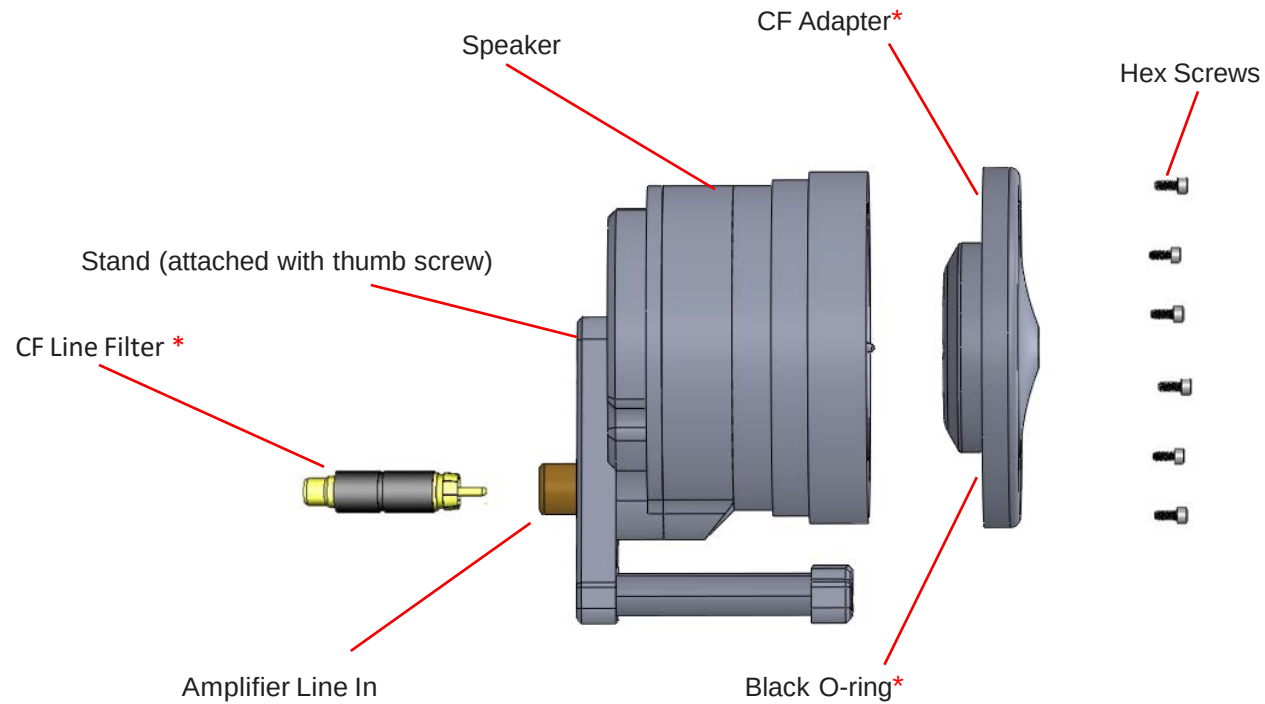

***Important: Closed Field Adapter, and CF Line Filter used only for Closed Field Configuration!**

Free Field Configuration



Closed Field Configuration

